


# SureForm from da Vinci Stapling

Achieve more with your stapler.

At its core is SmartFire technology, which monitors tissue compression before and during firing. With more than 1,000 measurements per second,<sup>1</sup> SmartFire makes automatic adjustments to optimize the staple line and prevent tissue damage, across a range of tissue thicknesses.<sup>1</sup> A 120° cone of articulation gives you the freedom to approach tissue where it naturally lies, potentially reducing the need to pull it towards the stapler.<sup>1</sup> When you choose SureForm staplers, you can potentially reduce uncertainty and increase stapling performance.<sup>1</sup>

SureForm staplers are available in 45 mm and 60 mm staple line lengths. SureForm 45 mm and 60 mm reloads include six rows of staples and a new knife with every reload.

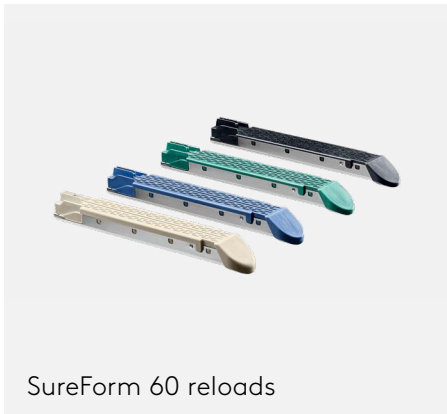


SureForm 60 reloads

- White, Blue, Green and Black reloads
- 6 rows of staples
- New knife with every reload



SureForm 60

- 60 mm staple line / 55 mm cut line
- An anvil centric design
- Straight-tip anvil
- Operates through a 12 mm cannula
- 12 fires per device
- Single patient use

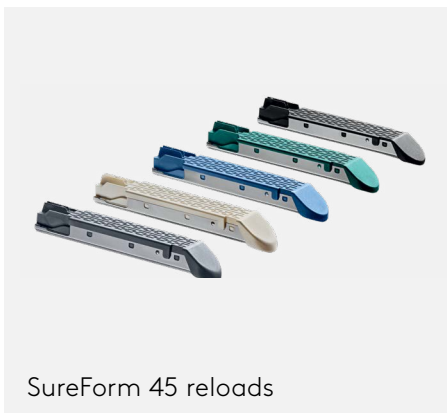

SureForm 45 reloads

- Gray, White, Blue, Green and Black reloads
- 6 rows of staples
- New knife with every reload



SureForm 45

SureForm 45 Curved-Tip

- 45 mm staple line / 40 mm cut line
- An anvil centric design
- Straight-tip or curved-tip anvil
- Operates through a 12 mm cannula
- 12 fires per device
- Single patient use

1. Results from internal testing on file.

---

## Ordering Information

### SureForm stapler

---

Part Number	Product Description	Uses	Quantity/Box
480460	SureForm 60 stapler	Single	6
480445	SureForm 45 stapler	Single	6
480545	SureForm 45 curved-tip stapler	Single	6

### SureForm stapler reloads

---

Part Number	Product Description	Uses	Quantity/Box
48360W	SureForm 60 white reload (2.5 mm, 6-row)	1 firing	12
48360B	SureForm 60 blue reload (3.5 mm, 6-row)	1 firing	12
48360G	SureForm 60 green reload (4.3 mm, 6-row)	1 firing	12
48360T	SureForm 60 black reload (4.6 mm, 6-row)	1 firing	12
48345M	SureForm 45 grey reload (2.0 mm, 6-row)	1 firing	12
48345W	SureForm 45 white reload (2.5 mm, 6-row)	1 firing	12
48345B	SureForm 45 blue reload (3.5 mm, 6-row)	1 firing	12
48345G	SureForm 45 green reload (4.3 mm, 6-row)	1 firing	12
48345T	SureForm 45 black reload (4.6 mm, 6-row)	1 firing	12

---

## Ordering Information (continued)

### 12mm stapling accessories

Part Number	Product Description	Uses	Quantity/Box
470375	12 mm and stapler cannula (100 mm)	NA	1
470376	12 mm and stapler blunt obturator	NA	1
470389	12 mm and stapler cannula, long (150 mm)	NA	1
470390	12 mm and stapler blunt obturator, long	NA	1
470380	12 mm and stapler cannula seal (box of 10)	1 procedure	10
470381	12 – 8 mm reducer (box of 6)	1 procedure	6
470395	12 mm and stapler bladeless obturator	NA	1
470396	12 mm and stapler bladeless obturator, long	NA	1

#### To learn more about da Vinci Stapling products, visit:

[www.intuitive.com/en-gb](http://www.intuitive.com/en-gb)

#### To contact a representative or receive additional information:

Visit [www.intuitive.com/en-gb](http://www.intuitive.com/en-gb) or call Intuitive Customer Service:

#### Europe

+41.21.821.2020 or  
+800.0821.2020 (toll free)

#### Product Information

##### SureForm 45 and 60

The SureForm 60 and SureForm 45, SureForm 60 and SureForm 45 Reloads are intended to be used with a compatible da Vinci Surgical System for resection, transection, and, or creation of anastomoses in General, Thoracic, Gynecologic, Urologic, and Pediatric surgery. The devices can be used with staple line or tissue buttressing material (natural or synthetic).

The SureForm 60 Stapler and SureForm 45 Stapler, SureForm 60 and SureForm 45 Stapler Reloads are class IIa and IIb medical devices CE marked (CE 2460) under the European Medical Devices Directive (93/42/EEC), manufactured by Intuitive Surgical, Inc. Refer to Instructions For Use before use.

##### Legal Notices

Individuals' outcomes may depend on a number of factors, including but not limited to patient characteristics, disease characteristics and/or physician/surgeon experience.

The information contained in this document has been checked and compiled with the greatest care. However, no responsibility is taken for its correctness, completeness and topicality. It is the sole responsibility of the recipient to check all information before using it in the individual case.

Some products, features or technologies may not be available in all countries. Please contact your local Intuitive representative for product availability in your region.

Refer to the product specific User Manual for indications, contraindications, warnings and other product information.

##### Privacy Notice

Intuitive's Privacy Notice is available at [www.intuitive.com/privacy](http://www.intuitive.com/privacy).

© 2022 Intuitive Surgical Operations, Inc. All rights reserved. Product and brand names/logos are trademarks or registered trademarks of Intuitive Surgical or their respective owner.